

Rec'd PCT/PTO

24 JAN 2005

REC'D 17 OCT 2003  
WIPO PCT

PI 1076314

# THE UNITED STATES OF AMERICA

TO ALL TO WHOM THESE PRESENTS SHALL COME:

UNITED STATES DEPARTMENT OF COMMERCE  
United States Patent and Trademark Office


October 09, 2003

THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM  
THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK  
OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT  
APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A  
FILING DATE.

APPLICATION NUMBER: 60/400,332

FILING DATE: August 01, 2002

RELATED PCT APPLICATION NUMBER: PCT/US03/23928

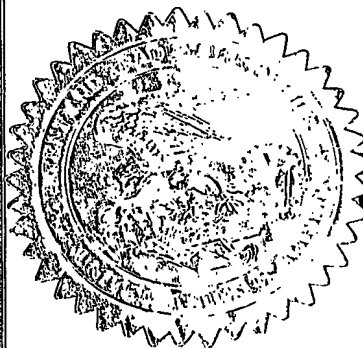

By Authority of the  
COMMISSIONER OF PATENTS AND TRADEMARKS

*M. Tarver*

M. TARVER  
Certifying Officer

**PRIORITY  
DOCUMENT**

SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH RULE 17.1(a) OR (b)

BEST AVAILABLE COPY

08-02-01

08/01/02  
 SUBSTITUTE for Provisional Application for Patent Cover Sheet PTO/SB/16 (10-01)  
 Approved for use through 10/31/2002 OMB 0651-0032  
 U.S. Patent and Trademark Office, U.S. DEPARTMENT OF COMMERCE

# PROVISIONAL APPLICATION FOR PATENT COVER SHEET

This is a request for filing a PROVISIONAL APPLICATION FOR PATENT under 37 CFR 1.53 (c).

		DOCKET NUMBER		20930PV	
INVENTOR(S)					
Given Name (first and middle [if any])		Family Name or Surname		Residence (City and either State or Foreign Country)	
Michael J. Xu Su		Caulfield Liu Wang		Fort Washington, Pennsylvania 19034 North Wales, Pennsylvania 19454 Schwensville, Pennsylvania 19473	
<input type="checkbox"/> Additional inventors are being named on the separately numbered sheets attached hereto					
TITLE OF THE INVENTION (500 characters max)					
NOVEL PROCEDURE FOR GROWTH, IMAGING, AND ENUMERATION OF MICROBIAL COLONIES FOR SEROLOGICAL OR SCREENING ASSAYS					
CORRESPONDENCE ADDRESS					
Direct all Correspondence to: Merck & Co., Inc. Patent Department - RY60-30 P.O. Box 2000 Rahway					
STATE		New Jersey		ZIP CODE	07065
COUNTRY		U.S.A.			
ENCLOSED APPLICATION PARTS (check all that apply)					
<input checked="" type="checkbox"/> Specification		Number of Pages		36	
<input type="checkbox"/> Drawing(s)		Number of Sheets		10	
<input type="checkbox"/> Application Data Sheet. See 37 CFR 1.76					
<input type="checkbox"/> CD(s), Number					
<input type="checkbox"/> Other (specify)					
METHOD OF PAYMENT OF FILING FEES FOR THIS PROVISIONAL APPLICATION FOR PATENT (check one)					
<input type="checkbox"/> A check or money order is enclosed to cover the filing fees					
<input checked="" type="checkbox"/> The Commissioner is hereby authorized to charge filing fees or credit any overpayment to Deposit Account Number:				13-2755	
FILING FEE AMOUNT (\$)				\$160.00	

The invention was made by an agency of the United States Government or under a contract with an agency of the United States Government.

☐ No.

☐ Yes, the name of the U.S. Government agency and the Government contract number are: \_\_\_\_\_

Respectfully submitted,

SIGNATURE \_\_\_\_\_

Date 08/01/2002

TYPED or PRINTED NAME Anna L. Cocuzzo

REGISTRATION NO 42,452  
(if appropriate)

TELEPHONE 732-594-1273

## EXPRESS MAIL CERTIFICATE

DATE OF DEPOSIT August 1, 2002

EXPRESS MAIL NO. EL523911511US

I HEREBY CERTIFY THAT THIS CORRESPONDENCE IS BEING DEPOSITED WITH THE UNITED STATES POSTAL SERVICE AS EXPRESS MAIL "POST OFFICE TO ADDRESSEE" ON THE ABOVE DATE IN AN ENVELOPE ADDRESSED TO ASSISTANT COMMISSIONER FOR PATENTS, WASHINGTON, D.C. 20231

MAILED BY

Christina Caffrey

DATE

8-1-02

In Duplicate

## TITLE OF THE INVENTION

NOVEL PROCEDURE FOR GROWTH, IMAGING, AND ENUMERATION OF  
MICROBIAL COLONIES FOR SEROLOGICAL OR SCREENING ASSAYS

## 5 BACKGROUND OF THE INVENTION

Microorganisms and the impact they exert on resident hosts are of great concern and importance and the continued study thereof commands great resources. Microorganisms, while indispensable components of our ecosystem and significant contributors to the production of various important consumer products such as antibiotics, vaccines, vitamins, and enzymes, also have very negatively impacted humans and disrupted society over the millennia. For this reason, one essential area of focus in the area of microbiology is in the evaluation of various treatment regimens (e.g., antibiotics) on the viability of various microorganisms; e.g., bacteria, viruses, algae, fungi, and protozoa.

15 Bacterial pathogens, for instance, are continuously evolving a more refined and increasingly powerful resistance towards proven-effective antibacterial drugs. This trend and the resulting drug- and multi-drug-resistant pathogens have become more of an issue in the last decade. As a result, efforts have largely focused on identifying alternatives to the current arsenal of available antibacterial drugs.

20 Vaccination through the elicitation of a host immune response offers a very effective alternative to administering antibacterial drugs and has remained the focus of ongoing antibacterial efforts. Antibacterial vaccines are being developed at a very rapid rate. More, assays useful in evaluating the efficacy of developed vaccines and various antimicrobial agents are continuously being scrutinized and further optimized. This latter focus of optimization of assays and, in a preferred embodiment, antibacterial assays is advanced herein. Accordingly, a brief discussion of two basic types of assays utilized in the evaluation of antibacterial agents follows.

25 The first type of assay detects specific antibodies generated in response to a bacterial antigen or an administered vaccine. Assays capable of functioning in this capacity include radioimmunoassays (RIAs) and enzyme-linked immunosorbent assays (ELISAs); see, e.g., Schiffman *et al.*, 1980 *J. Immunol. Methods* 33:133-44; Nahm *et al.*, 1996 *J. Infect. Dis.* 173:113-118; Quataert *et al.*, 1995 *Clin. Diagn. Lab. Immunol.* 2:590-597. Briefly, assays of this nature are directed at measuring total binding and do not provide information regarding functionality of the binding antibodies.

35

The proper functioning of specific antibodies in question is investigated by a second type of assay capable of measuring the efficacy of the antibodies to bind a specific antigen or agent and mark the antigen/agent for eventual destruction by complement in the absence or presence of phagocytic cells. These  
5 assays play a very important role in predicting antibiotic efficacy.

An assay in accordance with this description is the opsonophagocytic assay. Opsonophagocytosis is, briefly, a process whereby an invading cell or microbe, e.g., *Streptococcus pneumoniae*, is bound by circulating antibody and complement (or complement components). The bound antibody activates the complement cascade  
10 resulting in the deposition of certain complement components (e.g. C3b) onto the surface of the microbe. Phagocytic cells or alternative effector cells that have Fc receptors or complement receptors can then engulf and kill the opsonized microbe. Opsonophagocytic assays which are designed to mimic this process have traditionally employed peripheral blood leukocytes (PBLs) as effector cells and generally measure  
15 activity by a variety of techniques including radioisotopic, flow cytometric, microscopic, and viability assays; see, e.g., Obaro *et al.*, 1996 *Immunol. Lett.* 49:83-89; Vioarsson *et al.*, 1994 *J. Infect. Dis.* 170:592-599; Lortan *et al.*, 1993 *Clin. Exp. Immunol.* 91:54-57; Kaniuk *et al.*, 1992 *Scand. J. Immunol.* 36 (Supp. 11):96-98; Esposito *et al.*, 1990 *APMIS* 98:111-121; Sveum *et al.*, 1986 *J. Immunol. Methods*  
20 90:257-64; Guckian *et al.*, 1980 *J. Infect. Dis.* 142:175-190; and Winkelstein *et al.*, 1975 *Proc. Soc. Exp. Biol. Med.* 149:397-401. A standardized form of this assay was demonstrated wherein culturable phagocytes (differentiated HL-60 cells) were used to measure complement-dependent opsonophagocytic activity in sera from individuals vaccinated with various pneumococcal vaccines; Romero-Steiner *et al.*, 1997 *Clin.*  
25 *Diagn. Lab. Immunol.* 4(4):415-422. The use of culturable cells (like those of Romero-Steiner) was disclosed to eliminate the need for human donors and decrease the interassay variability that occurs with random PBL donors. *Id.*

There are alternative means for screening antibacterial agents for effectiveness against the growth and/or viability of a particular microorganism. The  
30 serum bactericidal assay (SBA) is one other example of a functionality-based assay capable of evaluating antibodies produced in response to an administered vaccine or bacterial antigen of interest. The purpose of the SBA assay is basically to evaluate whether anticapsular antibodies produced to a bacterial antigen, in combination with complement, are sufficient to confer host protection against invasive disease. This  
35 assay has proven particularly effective in the evaluation of vaccines directed against

*Neisseria meningitidis* infection, wherein circulating antibody and complement have been shown to confer host protection to meningococcal disease as early as 1918; Kolmer *et al.*, 1918 *J. Immunol.* 3:157-175. SBA titer is currently used as a standard by which to determine the effectiveness of proposed vaccine candidates in phase I and II field trials. SBA titer also serves to indicate seroconversion after immunization with currently licensed polysaccharide vaccines; *see, e.g.*, Anderson *et al.*, 1994 *Infect. Immunol.* 62:3391-3395; Milagres *et al.*, 1994 *Infect. Immunol.* 62:4419-4424; Zangwill *et al.*, 1994 *J. Infect. Dis.* 169:847-852. Accordingly, the assay is highly valued in the vaccine and antibiotic area. A standardized SBA assay has been developed at the Center for Disease Control; *see* Maslanka *et al.* 1997 *Clin. Diagn. Lab. Immunol.* 4:156.

Microbial assays such as the bacterial assays mentioned above serve a fundamental purpose in the research evaluation of antimicrobial agents and form a major part of any antimicrobial clinical program. For this reason, it is of great importance to continually design more efficient and optimized means by which to run these and other similar assays. In pursuit of this goal, Applicants have identified the following area in these various types of assays as ripe for improvement.

As one of skill in the art will appreciate, the above assays and microbial assays in general are traditionally run on multi-well plates and generally involve the transfer, and spreading, of an aliquot from each serum dilution of a treated sample to a forum for growth (typically, an agar petri plate). The bacteria or tested microorganism(s) are allowed to colonize on the plate and the appropriate analyses are subsequently carried out so that a determination can be made as to the effectiveness of the purported antimicrobial agent on the growth and/or viability of the microbe. One noted improvement in the art was from spreading aliquots of bacteria, for example, to a more refined "spotting" of bacterial samples (e.g., 5 microliter samples) onto small sections of an agar plate. This allowed for the accommodation of up to 48 samples on a single agar plate.

A common element in these assays, however, remains a limiting factor. Colonies are either counted manually with the aid of a microscope or counted via automated means *per agar plate*; processes, no doubt, burdened with limited throughput, very labor-intensive, timely, and fraught with opportunities for error. When it is considered that a clinical trial testing a multivalent vaccine may require performing several thousand functional antibody assays, this is a very significant limitation.

A more efficient method for performing the required assays and documenting the results would clearly be advantageous.

#### SUMMARY OF THE INVENTION

5           The instant invention relates to a novel method, for use within serological or screening assays, wherein microbes are grown as colonies on filter membranes in multi-well plates, according to the following process. A sample containing a given microbe (bacteria, for instance) in a liquid (or other transferable) medium, is added to the wells of a multi-well filter plate, for example a Millipore™  
10 Multiscreen™ 96 well filter plate. Excess medium is then removed by a process of vacuum filtration, centrifugation or other suitable means and a nutrient source (in the form of a growth medium) is provided (e.g., THYE broth). Importantly, residual growth medium trapped in or under the filter membrane enables the growth of microbes in discrete colonies on the surface of the filter. Growth of the microbes in  
15 this manner allows for the colonies to be stained, imaged, and counted automatically using such automated systems as, e.g., computer and video-based imaging systems. The assay can, thus, be exploited in evaluating the effectiveness of various antimicrobial agents on the growth and/or viability of various microorganisms in a timely and efficient manner.

20

#### BRIEF DESCRIPTION OF THE DRAWINGS

FIGURE 1 illustrates opsonophagocytic assay analyses employing the methods described herein wherein the subject bacteria was *Streptococcus pneumoniae* serotype 14.

25

FIGURE 2 illustrates colony growth obtained upon employing opsonophagocytic assay analyses according to the methods described herein wherein the subject bacteria was *Streptococcus pneumoniae* serotype 14.

FIGURE 3 illustrates colony counts obtained upon employing opsonophagocytic assay analyses according to the methods described herein wherein  
30 the subject bacteria was *Streptococcus pneumoniae* serotype 14.

FIGURE 4 illustrates serum bactericidal assay analyses employing the methods described herein wherein the subject bacteria was *Neisseria meningitidis* serotype C.

FIGURE 6 illustrates, in tabular format, results obtained upon serum  
5 bactericidal assay analyses according to the methods described herein wherein the  
subject bacteria was *Neisseria meningitidis* serotype C.

FIGURE 8 illustrates the experimental layout for the specificity studies of the analytical validation experiments. Specificity was assessed by determining the ability of polysaccharides of a known serotype (6B, 9V, 14, 18C, 19F, 23F, and C-Ps at 1 µg/ml) to inhibit killing by positive control antisera (Pool 019 and QC-1 tested at the 1:64 dilution).

FIGURE 9 illustrates the results obtained upon utilizing the methods disclosed herein in a serum bactericidal assay to screen hybridoma culture supernatants for the selection of functional monoclonal antibodies that kill *Neisseria meningitidis* serotype B.

The instant invention relates to a novel method wherein microbes such as bacteria are grown on filter membranes in multi-well plates, enabling the ready analysis of colonies grown thereon. Through the use of this process, the colonies can be readily fixed (i.e., killed), stained, imaged, and counted automatically using such automated systems as, *e.g.*, computer and video-based imaging systems. As such, the instant process avoids the manual counting of colonies grown on agar plates, timely and labor-intensive analyses, and the opportunity for error presented thereby. Also,

since microbes in the colonies are killed (fixed) in the process of staining, the method has safety advantages when working with pathogenic microorganisms.

Intrinsic to the instant discovery is the fact that Applicants had realized that residual growth medium trapped in or under a filter membrane is sufficient to enable microbes on the surface of the filter to grow as discrete colonies. The growth obtained on the filter plates is, further, consistent with that seen on agar for purposes of screening downstream antimicrobial activity and, importantly, does not detract from assay precision.

This was not appreciated prior to the instant invention. Filter membranes, to Applicants' knowledge, have been employed in the capture and isolation of microbes, and generally as a stage for various screening and serological assays, such as those carried out in multi-well plates. Upon completion of the reaction(s), however, the bacteria or tested microbe were transferred to a growth medium (*e.g.*, an agar petri plate), and manually counted in order to derive any experimental conclusions. Although devices exist for counting colonies on large (*e.g.*, 100 mm diameter) petri dishes, these devices are not generally amenable to rapid and automated counting of colonies grown in a multi-well (*e.g.*, 48 or 96 well) format. Efficient imaging and counting systems for assays run in multi-well plate format (*e.g.*, the C.T.L. imaging system for ELISPOT assay) have relatively recently become available. These imaging technologies, however, are not capable of enumerating microbial colonies grown on agar due to the irregularities of the agar surface and the difficulties in uniformly plating multiple samples of microbes (*e.g.*, bacteria) on a single plate.

The instant invention bypasses the art in that it provides a means of employing the filter membrane as a forum for microbial growth in a format amenable to automated analyses and avoids, in part, the agar petri plate which does not currently lend itself to efficient automated analyses in a multi-well format. Automated counting of samples in multi-well format has proved more consistent than the manual method, and the assay has proven acceptably rugged to changes in cell passage, operator, plate and counting method; Example 5. There is, further, no evidence that the counting method within the plates affected titer. Notably, individual results tended to be more variable on the agar plate than on the well plate.

The method is, therefore, particularly useful for various serological or screening assays in which determining the effect (microbicidal or microbistatic) of a biological or chemical agent on the number of a particular microorganism of interest



microbes is desired. In certain applications, the reaction between the test agent and the microbes can take place in a separate reaction container. The process then generally involves (1) transferring a sample comprising the microbe (*e.g.*, bacteria) in a liquid (or other transferable) medium to the wells of a multi-well filter plate (*e.g.*, a Millipore™ 96 well filter plate or an equivalent thereof (an equivalent thereof being defined as a plate in a multi-well format comprising a filter compatible for use therein)); (2) removing excess media (media other than that captured within and/or under the filter membrane) by a process of vacuum filtration, centrifugation or other means suitable for removing liquid (or other) medium from multi-well plates; and (3) allowing sufficient time for the microbes (*e.g.*, bacteria) to grow into discrete colonies for subsequent analyses. In alternative applications, the assay could feasibly be carried out entirely within and on the microbial filter plate.

Filter plates of use in the instant invention are those suited for use in a multi-well format. The term "filter plates" employed throughout the instant application is to be interpreted as including both specifically crafted "filter plates" in multi-well format as well as simply 96 well plates comprising filters. Particularly preferred are Millipore™ 96 well plates, *e.g.*, the 0.45µm (pore size) Durapore PVDF filters. Most preferred embodiments of the instant invention employ the Millipore™ MultiScreen™ 96 well plates. It is to be noted that the instant invention is not limited to wells contained within a 96 well format. Any multi-well format suited for ready analysis via automated means is definitely encompassed hereby. Such capabilities are enabled by Applicants' finding that microbes can be manipulated to grow on filter membranes in multi-well format in distinct colonies to an extent comparable for vaccine evaluation purposes to microbes grown on agar.

Filter plates possessing the following characteristics (alone or in combination) have been decidedly preferred for use in the disclosed methods: low levels of protein binding, compatibility with bacterial growth, sterility, and hydrophilicity. Plates possessing all four of these characteristics are most preferred for use within the instant invention. Particularly preferred are Millipore™ HV plates. Opaque plates are further preferred for use within the instant invention as they retard light refraction when undergoing automated analyses. Specifically preferred embodiments of the instant invention employ Millipore™ MultiScreen™ HV 0.45 µm Opaque Sterile Filtration Plates.

Filter plates provide a forum for growth within the well. Bacteria remain on the surface of the filter plates as they are larger than the pore size of the

filter membrane. Excess media retained within the multi-well plates is then removed by either vacuum filtration, centrifugation or other means found suitable for removing liquid or alternative transferable medium from multi-well plates. The medium provided to the microbial sample following transfer to the multi-well plate can be any  
5 nutrient medium (growth medium) provided to the bacteria or tested microbe(s) following transfer to the filter plates.

The filters, despite the removal process, will retain some medium within or under the filter membrane. Applicants have discovered that this residual medium is surprisingly sufficient to provide adequate nutrients to enable and support  
10 the growth of microbes into discrete colonies on the surface of the membrane. Further and importantly, this growth is sufficient for purposes of evaluating the effectiveness of antimicrobial (and, preferably, antibacterial) vaccine candidates. A comparison of the results obtained upon employing bacteria grown on an agar medium with the results obtained upon the use of bacteria grown on filter membranes in the wells of a  
15 multi-well plate showed no demonstrable impact of the particular growth forum. Quite to the contrary, analysis in 96 well format proved more consistent, the assay proved acceptably rugged to changes in cell passage, operator, plate and counting method, and there was no evidence that the counting method within the plates affected bacterial titer; Example 5.

20 Removal of the medium is noted to be important to the disclosed methods. Applicants have found that medium beyond residual medium (that trapped within or under the filter) is undesirable. Insufficient removal results in growth of the microbe(s) as a homogeneous suspension, rather than as discrete colonies required for accurate enumeration and analysis. More, it is preferred that the medium is removed  
25 in such a manner as to increase the contact of microbes present in the sample with the filter membrane (*e.g.*, as in vacuum filtration and centrifugation, where the liquid is funneled downwardly through the membrane).

As indicated above, following removal of medium from the wells, the microbe(s) are incubated for an amount of time sufficient to permit growth into  
30 discrete colonies for subsequent analyses; preferably, 14-18 hours, but generally dependent on the particular microorganism. Preferably, the plates comprising the microbes are kept hydrated. This can be accomplished through a number of means such as incubation in a humidified environment such as a water-jacketed humidified incubator or by covering the filter plate in, for instance, a Ziploc bag. Any suitable  
35 alternative serving to accomplish this same function is also encompassed hereby.

An interesting feature of the instant method is that the colony size is limited by the available nutrients trapped within or under the filter membrane. Thus, when there are few colonies, the individual colonies are relatively large, whereas when there are many colonies, the individual colonies are smaller. The benefit of this particular finding is that, in certain applications, up to approximately 300 colonies per well can be enumerated in an individual well, rendering the dynamic range of the method wider than that for the previously employed agar-based methods in which the maximum number of colonies discernable after "spotting" of 5 microliter samples is approximately 70.

After some time during which the microbe(s) are allowed to grow into discrete colonies on the filter plate, the colonies are fixed and stained; a preferred stain of which is Coomassie blue, but which can be any agent capable of providing a means for specifically highlighting an object for detection by the visualization, imaging and/or enumeration system employed. Staining in general terms enables ready visualization of the colonies by the system employed. The colonies are then analyzed. This is accomplished with any device suited for analysis of 96 well plates, or whichever multi-well format is utilized. Preferred for use in the instant invention is any computer-assisted video imaging and analysis system. In a particularly preferred embodiment, the computer-assisted imaging and analysis system is the ImmunoSpot™ Analyzer offered by C.T.L. (Cleveland, OH), or a similar imaging system offered by, for instance, Resolution Technology, Inc. (Columbus, OH).

Applicants believe the methods disclosed herein are widely applicable to various microorganisms, e.g., bacteria, viruses, algae, fungi (e.g., *Candida albicans* and *Aspergillus*) and protozoa. Samples can be obtained from any of a number of sources. For purposes of exemplification, the bacteria can be isolates of *Streptococcus pneumoniae*, *Neisseria meningitidis*, *E. coli*, *Staphylococcus aureus*, *Bacillus anthracis*, or any gram-positive or gram-negative bacteria.

The bacteria or other microbe(s) are generally contained in a liquid (or other transferable) medium. A medium containing a nutrient source (a growth or nutrient medium) must be provided. Broth as a nutrient medium is particularly preferred. Particular embodiments of the instant invention employ Todd-Hewitt ("TH") yeast extract ("THYE") broth. In preferred embodiments, said broth is employed for the growth of Streptococci. Alternative embodiments employ tryptic soy broth ("TSB"). In specific embodiments, said TSB broth is utilized for the growth of *Neisseria* and *E. coli*. Alternative broth media suitable for the growth of the

microorganism of interest can be used. In the examples provided, 100  $\mu$ l per well was added to the Millipore™ 96 well HV plate.

Using the disclosed methods, one of skill in the art can evaluate the efficacy of any particular antimicrobial agent. The antimicrobial agent of interest which can be any compound (*e.g.*, antibiotic or microbistatic agent) or biological (*e.g.*, antiserum) having an impact on the growth and/or viability of the microbe(s) of interest is placed in contact with the microbial sample and left for a period of time. Transfer of the microbial sample to the multi-well plate can take place before or following contact with the antimicrobial agent. In the situation wherein the microbial sample is brought into contact with the antimicrobial agent prior to transfer, this contact takes place in a different medium, for instance another well or an agar plate. In specific embodiments, the microbial samples to be analyzed are transferred from mediums wherein various functional assays, such as opsonophagocytic and serum bactericidal assays were run. In one specific embodiment, prior to transfer to the filter plate, the bacteria are "preopsonized" (*i.e.*, placed in contact with antibody (*see, e.g.*, Example 2D)); and placed in contact with complement (or the active components thereof) and differentiated HL-60 cells (or any cells capable of clearing the marked bacteria, *e.g.*, polymorphonuclear leukocytes). Antibody which can be utilized in these experiments can be derived from an individual or animal which was administered the microbe or an effective antigen of same. Antibodies specific to the bacteria or microbe of interest contained within the sample can also be obtained commercially or via generation outside of the human (*e.g.*, in rabbit). In another embodiment of the instant invention, prior to transfer to the filter plate, the bacteria are preopsonized and then placed in contact with complement (or active components thereof); (*see, e.g.*, Example 4D).

Antimicrobial agents which can be evaluated via the utilization of the methods disclosed herein can be selected from any compound (*e.g.*, antibiotic or microbistatic agent) or biological (*e.g.*, antiserum) that impacts the growth and/or viability of the specific microbe of interest. Hybridoma culture supernatant can be screened for the production of monoclonal antibodies effective against a microbe of interest.

The following non-limiting Examples are presented to better illustrate very specific embodiments of the present invention. As one of ordinary skill in the art

will appreciate, there are numerous ways in which one can carry out the various methods disclosed herein and accomplish the same or similar results.

#### EXAMPLES

- 5                   The required precautions and procedures outlined for compliance with OSHA's Blood Borne Pathogen Standard should be followed when handling human blood products.

#### 10       EXAMPLE 1. MATERIALS AND EQUIPMENT FOR THE OPSONOPHAGOCYTIC ASSAY

##### A. Materials

1. 96-Well Cell Cluster Plates (U-bottom) (Costar, cat no. 3799)
2. Reagent Reservoirs, (Fisher Scientific, cat no. 13-681-101)
3. Sterile (0-200  $\mu$ l) pipette tips
- 15 4. Tips for Matrix Pipette (Matrix Technologies, 0-850  $\mu$ l, cat no. 8052)
5. Serology disposable pipettes, sterile (2, 5, 10, 25, 50 mL volume)
6. Hemocytometer or equivalent (e.g. Kova Glasstic Slide (Fisher Scientific, cat no. 22-270141))
7. Centrifuge tubes: 15 ml and 50 ml conical polystyrene centrifuge tubes
- 20 8. Millipore™ MultiScreen™ filtration plate, HV, Opaque sterile Cat # MHVBS4510
9. Whatman filter paper No. 5, circles, 150mm diameter, Cat # 1005 150
10. CD-Recordable Disks
11. Ziploc Bags for HV plate overnight incubation.
12. HV plate storage, vacuum-sealable bags

##### 25       B. Equipment

1. CO<sub>2</sub> Incubator set at 37 $\pm$ 2 °C, 5 $\pm$ 1% CO<sub>2</sub>
2. Centrifuge for cell washing (Jouan CR422)
3. Matrix Pipette (Matrix Technologies, 15-850  $\mu$ l, cat no. 2014)
4. Single channel pipettes: 1-10 $\mu$ l, 20-100 or 200 $\mu$ l, and 100-1000 $\mu$ l
- 30 5. Multi channel pipettes (5-50  $\mu$ l, 50-200  $\mu$ l)
6. Inverted microscope with 1X and 10X objectives
7. Compound microscope
8. Horizontal rotator shaker (up to 220 rpm range) Thermolyne RotoMix type 50800
9. Colony counting system: ImmunoSpot™ counter (Cellular Technology Limited
- 35 (C.T.L.), Cleveland, OH).

C. Reagents

1. Gelatin (2 % solution) (Sigma, cat no. G1393)
2. Hank's Buffered Saline Solution (HBSS) (Gibco-BRL Life Technologies, cat no. 14025-092)
- 5 3. Hank's Buffered Saline Solution (HBSS) without  $\text{Ca}^{2+}$ ,  $\text{Mg}^{2+}$  (Gibco-BRL Life Technologies, cat no. 14175-095)
4. Todd-Hewitt Broth (Becton Dickinson, cat no. 8117360, Fisher Scientific, cat no. B11736)
5. Yeast Extract powder (#28926), Difco
- 10 6. 3-4 week old rabbit complement, sterile (Pel. Freez, code no. 31061)
7. HL-60 cells differentiated as described in Romero-Steiner *et al.*, 1997 *Clin. Diagn. Lab. Immunol.* 4(4):415-422.
8. 0.4% trypan blue solution in PBS, pH 7.2
9. *S. pneumoniae* bacteria (frozen aliquots of various serotypes: 1, 4, 6B, 9V, 14, 15 18C, 19F, and 23F)
10. *S. pneumoniae* capsular polysaccharides (various serotypes: 4, 6B, 9V, 14, 18C, 19F, and 23F)
11. 0.01% Coomassie blue prepared in 45% methanol, 45% water and 10% acetic acid
- 20 12. Coomassie blue diluent: 45% methanol, 45% water and 10% acetic acid
13. Sterile water (RCM-66)

D. Media

1. Todd-Hewitt Broth with 0.5% yeast extract, 50ml aliquot, autoclaved.
2. Gelatin (0.1%) in HBSS with  $\text{Ca}^{++}$ ,  $\text{Mg}^{++}$  ("Opsono buffer")
- 25

EXAMPLE 2. PROCEDURE FOR OPSONOPHAGOCYTIC ASSAYPlate layout for OPK assay (serotype 6B)

Sample 1		Sample 2		Sample 3		Sample 4		Sample 5		Controls	
1	2	1	2	1	2	1	2	1	2	1	2
1:8		1:8		1:8		1:8		1:8		Cells + C'	
1:16		1:16		1:16		1:16		1:16		(no serum)	
1:32		1:32		1:32		1:32		1:32		QC serum (2X titer)	
1:64		1:64		1:64		1:64		1:64		QC serum (1X titer)	
1:128		1:128		1:128		1:128		1:128		QC serum (0.5X titer)	
1:256		1:256		1:256		1:256		1:256		QC sera + 6B PS*	
1:512		1:512		1:512		1:512		1:512		QC sera + C-PS*	
1:1024		1:1024		1:1024		1:1024		1:1024		Medium alone	

### A. Test Sample Dilutions

Note: Serum samples are generally tested without pre-dilution. A minimum of 40  $\mu$ l is recommended needed to test a serum sample (in duplicate) for opsonic activity against a single serotype of pneumococci. Per specific experiment needs, sera may be heat-inactivated ( $56 \pm 2$  °C,  $30 \pm 5$  minutes) prior to testing, and repeat cycles of freezing and thawing should be avoided. Samples are run in duplicate, so A1 and A2 are the same serum, A3 and A4 from a second serum sample, and so on until wells A11 through A12 are reached. These last two columns are used for the complement ("C") controls and for Quality Control ("QC") tests of assay validity (see plate template above).

1. A sterile reagent reservoir was filled with Opsono buffer and a multi-channel pipette was used to add 10  $\mu$ l of Opsono buffer to all wells in columns 1-10 except Row A.
2. 20  $\mu$ l of undiluted serum was aliquotted into the appropriate wells in Row A according to the plate set-up template.
3. Two-fold serial dilutions of the samples were made using a multichannel pipetter adjusted to 10  $\mu$ l by first transferring 10 $\mu$ l from wells A1 to A10 in the first row to the corresponding wells in Row B; avoiding air bubbles.
4. The contents of the wells were carefully mixed by pipetting up and down ~5 times.
5. 10  $\mu$ l from the second row (Row B) was transferred to the third row (Row C) and step 4 was repeated. This process was continued on through Row H. The excess 10  $\mu$ l from Row H was discarded into the waste.
6. All wells now contained 10  $\mu$ l. Note: Once the remainder of the reagents are added, Row A will become a 1:8 dilution. Row B = 1:16, row C = 1:32, row D = 1:64, row E = 1:128, row F = 1:256, row G = 1:512 and row H = 1:1024. Columns 11-12 will contain pre-diluted QC test sera and controls; see below.
7. 10  $\mu$ L of opsono buffer was added to all wells in columns 1-10 to bring the total volume to 20  $\mu$ L per well. Note: Buffer was not added to the wells in column 11 and 12 at that time.

### B. Plate quality control

Each plate had its own controls for assay validity. These controls were set up in columns 11 and 12:

1. Complement control wells were (A11, A12, B11, B12): 20 $\mu$ l of Opsono buffer was added to make the volume equal to the volume of the test sample wells.

2. Positive QC sera controls: QC sera were prepared by pooling serum from adult humans who had previously been vaccinated with one or more doses of pneumococcal vaccine. The serum pools were pre-tested in the OPK assay to determine the endpoint titer against each serotype. The QC serum was included in each assay plate at three dilutions expected to bracket the pre-determined titer for a given serotype. QC serum in appropriate wells was diluted to approximate dilutions of 2X (C11, C12), 1X(D11, D12), and 0.5X(E11, E12), the predetermined endpoint titer for that lot of QC serum. 10 µL of Opsono buffer was added to these wells, for a final volume of 20 µL.
3. Specificity controls:
- Serotype specific control (F11, F12): These wells contained 10µl of known titer sera which approximately represents the 2X concentration of OPK titer and 10µl of serotype-specific polysaccharide at 8µg/ml.
- Serotype non-specific control (G11, G12): These wells contained 10µl of known titer sera which give the 2X concentration of OPK titer and 10µl of non-specific bacteria cell wall polysaccharide (C-Ps) at 8µg/ml.
4. Medium alone control (H11, H12): 30 µl of Opsono buffer was added to these wells (H11 and H12).
5. For an assay plate(s) to be valid, it was noted that (1) the positive control QC serum must be positive (i.e., ≥50% inhibition) at least at one of three dilutions tested (e.g. 1:320, 1:640, or 1:1260; C11, C12; D11, D12; AND E11, E12, respectively); (2) The specificity control wells (quality control ("QC") serum (F11, F12)+ homologous polysaccharide ("PS"); G11, G12) must be negative (i.e., <50% inhibition); and (3) the average number of colonies across the four "No Serum" (complement control) wells (A11, A12, B11, B12) must be ≥25.
- C. Harvesting differentiated HL-60 Cells
1. Flasks were gently swirled to distribute cells and the contents of each flask were then poured into 50 mL centrifuge tubes.
2. The cells were centrifuged at 160 X g for 10 minutes at room temperature (Note: this step was started before adding bacteria to the plate for preopsonization).
- Preopsonization Steps 1-4 (Example 2D) were performed before proceeding to the next step.
3. Supernatant was poured into sterile bottles for disposal and the remaining medium was discarded by inverting tube with a single motion and draining excess medium by keeping the tube inverted until the supernatant was removed.



Note: Differentiated HL-60 cell medium contains DMF and should be disposed as chemical waste.

4. The HL60 cell pellet was resuspended from every 100 mL of original culture that was used in 40-50 ml of Hanks buffer (room temperature) without  $\text{Ca}^{++}$ ,  $\text{Mg}^{++}$  and phenol red. Note: If 200 mL of original cell culture are used, use 80-100 mL of the Hank's Buffer without  $\text{Ca}^{++}$ ,  $\text{Mg}^{++}$  and phenol red.
5. An aliquot of cells was removed and viable cell count was determined by staining cells with Trypan blue dye and counting cells using a hemocytometer or equivalent tool.
- 10 6. Meanwhile, cells were centrifuged at 160 X g for 10 minutes at room temperature and supernatant was removed as in step 3 of Example 2C. Note: Step 1 of Example 2E was performed at this point to allow the baby rabbit serum complement to thaw.
7. The HL-60 cell pellet was resuspended in opsono buffer at  $1 \times 10^7$  cells per ml. These cells were then ready to use in the assay, using 4 ml per plate. They were kept
- 15 at room temperature until their use. Note: A slight (e.g. at least 1 mL) excess of cells was purposely present for the assay. Note: only plastic pipettes were used. Cells attach to glass. Hanks was gently pipetted up and down 3 to 4 times only, and the pellet was dispensed slowly.

#### D. Preopsonization

- 20 1. A bacteria frozen stock of the desired strain from the  $-70^{\circ}\text{C}$  freezer was allowed to thaw at room temperature, and used immediately; otherwise, the thawed bacteria need be kept on ice until ready to make dilution but no longer than 1 hour.
2. The contents of the vial were mixed by slowly pipetting about one-half the volume of medium in the vial up and down 3 to 5 times; avoiding the introduction of air into
- 25 the medium. Note: Do not vortex. An initial dilution was made of 1:10 as well as subsequent dilutions as needed to reach the desired final concentration. Dilutions were made with Opsono buffer (usually in 2-3 steps, e.g., 1:10, then, 1:100 and, finally, 1:40 to make a 1:4000 dilution) in order to adjust the number of bacteria to ~1,000 cfu/10  $\mu\text{l}$ . Note: (Frozen aliquots of *S. Pneumoniae* were pre-tested to
- 30 determine the appropriate dilution.)
- Note: Each plate requires ~1 ml of diluted bacteria. Dilution varies according to each frozen strain and lot number. Make sure you have a slight (e.g. at least 1 mL) excess of cells for the assay.
3. A multichannel pipette set to 10  $\mu\text{l}$  was used to add the diluted bacterial suspension
- 35 to each of the 96 wells in a microtiter plate, except for the "medium alone" wells, H11

and H12. A sterile reagent reservoir was used to hold the diluted bacterial suspension. Where multiple plates were being run, the bacteria dilution was mixed by repeated aspiration and expulsion with a pipette, between additions to successive plates.

Volume in each well at this point:

	<u>Titer test</u>	<u>C' control</u>	<u>Specificity</u>	<u>Med. alone</u>
5 Undiluted or diluted serum:	10 $\mu$ l	-	10 $\mu$ l	-
Opsono buffer:	10 $\mu$ l	20 $\mu$ l	-	30 $\mu$ l
Polysaccharide (specif.'y. control):-	-	-	10 $\mu$ l	-
Bacterial suspension:	10 $\mu$ l	10 $\mu$ l	10 $\mu$ l	-
10 Total:	30 $\mu$ l	30 $\mu$ l	30 $\mu$ l	30 $\mu$ l

4. The plates were covered with their lids and the microtiter plates were incubated on a horizontal rotator shaker at 200 rpm at  $37 \pm 2$  °C,  $5 \pm 1\%$  CO<sub>2</sub> for  $15 \pm 2$  min. This step allowed the pre-opsonization of bacteria with diluted sera.

15 E. Complement Addition

1. Frozen stocks of complement were stored at -70°C. The vials of baby rabbit serum were thawed immediately before use. Note: The complement should be at  $\sim 4^\circ\text{C}$  before addition to plates. One vial of baby rabbit serum that contains at least 1 mL volume, per assay plate was thawed; and one additional vial was thawed to have a slight excess, to cover multiple plate usage.

20 2. Following the preopsonization, 10  $\mu$ l of complement source (Baby Rabbit Serum) was added to all 96 wells using a multi-channel pipetter.

25 3. When adding the complement, a 1ml vial was enough for one microtiter plate. All the vials necessary for the number of plates being run were thawed and combined in a sterile reagent reservoir.

4. The vials were pooled in a sterile reagent reservoir and the complement was quickly added to all wells.

F. Phagocytosis with HL-60 Cells

30 1. 40  $\mu$ l of the washed resuspended cells (at  $1 \times 10^7$  ml) was added into each well of a microtiter plate for the opsonophagocytic assay.

Note: The total number of cells added to each well were and should be  $4 \times 10^5$  in 40  $\mu$ l. Each 96 well plate requires 4 ml. A sterile reagent reservoir was used to hold the cell suspension.

Volume in each well at this point:

	<u>Titer test</u>	<u>C' control</u>	<u>Specificity</u>	<u>Medium alone</u>
Preopsonization	30 $\mu$ l	30 $\mu$ l	30 $\mu$ l	30 $\mu$ l
Baby rabbit serum	10 $\mu$ l	10 $\mu$ l	10 $\mu$ l	10 $\mu$ l
5 Differentiated HL-60 cells	40 $\mu$ l	40 $\mu$ l	40 $\mu$ l	40 $\mu$ l
<hr/>				
Total volume	80 $\mu$ l	80 $\mu$ l	80 $\mu$ l	80 $\mu$ l

2. The lids were replaced and the microtiter plates were incubated on a horizontal rotator shaker at  $220 \pm 20$  rpm at  $37 \pm 2$  °C,  $5 \pm 1\%$  CO<sub>2</sub>, for a period of  $45 \pm 5$  minutes.

10 G. Viable Bacteria Enumeration

1. Sterile water (RCM-66) was added to Millipore™ 96 well HV plate, 200 $\mu$ l per well.

2. The contents in the wells from the opsono reaction were carefully mixed (4-6 times) and the desired volume of sample (e.g. 5 or 10  $\mu$ l) from each well was

15 transferred into the corresponding wells of the HV plates using a multichannel pipette; this was done column by column.

3. Filtration: The HV plate was placed on the Millipore™ MultiScreen™ filtration system and the pump was turned on. Note: Make sure meter is showing  $10 \pm 2$  in. Hg-vac. The plate contents were filtered through the membrane by pushing down the plate with the lid to the base of the unit until the plate was fully emptied and no liquid remained visible.

4. THYE broth (room temperature) was added to the Millipore™ 96 well HV plate; 100 $\mu$ l per well.

5. Filtration: Step 3 was repeated.

25 Note: Filtration can happen very fast (within 5 – 10 seconds). Therefore, it is recommended that one monitor the filtration carefully to prevent over-drying of the plate. Note that too much broth remaining in the wells will give high background and poor quality staining. Remove plate from the manifold immediately after most of the broth is drawn through the wells, but avoiding over-drying and over-filtration.

30 Carefully monitor wells during filtration to observe when well contents have been aspirated, as discernable by the disappearance of liquid from the outside edges of the wells. Remove plate immediately after this occurs.

6. After filtration, the base of the HV plate unit was pat dry on a paper towel.

7. The plate was wrapped with plastic film to seal all moisture inside of the plate.

8. The wrapped plate was inserted in a sealable, plastic Ziploc bag and incubated overnight at  $37 \pm 2$  °C,  $5 \pm 1\%$  CO<sub>2</sub>, for a period of 20-24 hours. The plates were placed in an upright position in the incubator.

#### H. Colony staining and Counting

- 5 1. 0.01% Coomassie blue was filtered with Whatman No. 5 filter paper to remove undissolved Coomassie blue precipitates. Approximately 5ml of stain was required for each HV plate. Only fresh-filtered stain was used for HV plates.
2. Colony staining: The HV plate was placed on the Millipore™ filtration system without pushing down (Note: the plate was off the vacuum). The pump was turned  
10 on. Note: Make sure meter is showing  $10 \pm 2$  in. Hg-vac. 50µl per well of filtered Coomassie blue was added to the plate using Matrix Multichannel pipette. The plate was only stained for 15 to 30 seconds; then, the plate was pushed down onto the filtration manifold to empty the contents by vacuum filtration.
3. Destaining: Once the plate was fully emptied, 100µl per well Coomassie blue  
15 diluent was immediately added to each well and the plate was kept under the vacuum until no liquid was visible. This step removed most non-specific staining.
4. The base of the plate was pat dry on a paper towel or equivalent and the base was separated from the 96 well multiscreen plate.
5. The plate was dried upside down in a Biosafety hood with circulating air, until  
20 wells were completely dried; for approximately 10 minutes.
6. ImmunoSpot™ counter setting for colony counting

		Standard setting
	Sensitivity:	180
	Spot size:	0.03
25	Separation tolerance:	8.0
	Diffuse counting is:	on
	Diffuseness:	20
	Overdeveloped area handling is:	_____ off
	Color system:	Blue
30	Detailed counting is:	off
	Background balance is:	_____ on
	Background balance:	110
	Hole filling is:	on

Note: ImmunoSpot™ settings were calibrated and adjusted at 3 month intervals by  
35 comparison with counts obtained by human readers.

7. The plate was scanned in the imaging system, the image was printed and the data was counted. The image data were saved as tiff files on CD-R disk.

### I. OPK titer calculations

The phagocytic titer is the reciprocal of the serum dilution with at least 50% killing, when compared to the average growth in the complement control wells. Occasionally, sera with high titers need to be re-tested at higher initial dilutions than 1:8 to

- 5 determine the phagocytic titer. Serum samples with phagocytic titers <8 are reported as a titer of 4 for purposes of data analysis.

$$\% \text{ Killing} = \frac{\text{No Serum Control Avg CFU} - \text{Test Sample Avg CFU}}{\text{No Serum Control Avg CFU} - \text{Medium Only Avg CFU}} \times 100\%$$

### J. Experimental Results

- OPK assays were individually run with bacterial samples of *S. pneumoniae* serotypes  
10 6B, 9V, 14, 18C, 19F and 23F. Bacterial growth or lack thereof (due to the effect of an antibiotic substance) was able to be detected and analyzed within the multi-well plates. Figure 1 shows the basic assay template set-up utilized to run the experiment with serum samples of *S. pneumoniae* serotype 14. Figures 2 and 3, respectively, show the colony read-out obtained from the wells both pictorially and numerically.

15

## EXAMPLE 3. MATERIALS AND EQUIPMENT FOR SERUM BACTERICIDAL ASSAY

### A. Materials

1. 96-Well Cell Cluster Plates (Flat-bottom cell culture cluster, polystyrene) (Costar, cat no. 3595)
2. Reagent Reservoirs, (Fisher Scientific, cat no. 13-681-101)
3. Sterile (0-200 µl) pipette tips (Matrix)
4. Tips for Matrix Pipette (Matrix Technologies, 0-850 µl, cat no. 8052)
5. 14 ml polypropylene round-bottom tube, 17x 100 mm style (Becton-Dickinson, Falcon cat no. 352059)
- 25 6. 5 ml polystyrene round-bottom tube, 12x 75 mm style (Becton-Dickinson, Falcon cat no. 352058)
7. Millipore™ MultiScreen™ filtration plate, HV, Opaque sterile 0.45 µm hydrophilic, Durapore membrane. Cat # MHVBS4510
- 30 8. Whatman filter paper No. 5, circles, 150mm diameter, Cat # 1005 150
9. CD-Recordable Disks, Imation (ISG catalog # IMN-41181)
10. Fisherbrand Autoclave Bags for HV plate overnight incubation (8.5'' x 11'', Fisher Scientific #01-8151).

**B. Equipment**

1. CO<sub>2</sub> Incubator set at 37±2 °C, 5±1% CO<sub>2</sub>
2. Matrix Pipette (Matrix Technologies, 15-850 µl, cat no. 2014)
3. Single channel pipettes: 1-10µl, 20-100 or 200µl, and 100-1000µl (Lab Systems)
- 5 4. Multi channel pipettes (5-50 µl, 50-200 µl) (Lab Systems)
5. Horizontal rotator shaker (up to 220 rpm range) Thermolyne RotoMix type 50800
6. Colony counting system: ImmunoSpot™ counter (C.T.L., Cleveland, OH).
7. Millipore MultiScreen Separation System (Millipore, cat no. MAVM 096 0R)
8. Vacuum pump unit: VakuumSystem type ME 2S1 (Vacuubrand GMGH + Co)

**10 C. Reagents**

1. Glucose (Aldrich, cat no. 25307-3)
2. Dulbecco's PBS (D-PBS) (Sigma, cat no. D8537)
3. Tryptic Soy Broth (TSB): (Difco, Becton Dickinson cat no. 290612), sterile, 100 ml/bottle
- 15 4. 3-4 week old rabbit complement, sterile (Pel. Freez, code no. 31061)
5. Neisseria. Meningitidis bacteria (frozen aliquots of various serotypes: A, B, C, Y, and W-135): Originally obtained from Dr. Sandra Steiner, Centers for Disease Control and Prevention.
6. Neisseria. Meningitidis capsular polysaccharides (various serotypes: C, Y, and W-135) 1mg/ml in sterile water. Stored at -80 C. Merck Research Labs.
- 20 7. 0.01% Coomassie blue prepared in 45% methanol, 45% water and 10% acetic acid. Stored at room temperature for up to 3 months.
8. Sterile water (RCM-66)

**D. Media**

- 25 1. 10%glucose: prepared in sterile distilled water (RCM-66), 0.45µm filtered.
2. D-PBS-0.1% glucose (D-PBS-Glc): Dulbecco's PBS, 0.1% glucose, 0.45µm filtered. Prepared by 1:100 dilution of 10% glucose in D-PBS.

**EXAMPLE 4. PROCEDURE FOR SERUM BACTERICIDAL ASSAY**

Plate layout for SBA assay (serotype C)

	Serogroup C bacteria											
	Sample 1				Sample 2				Sample 3			
	Pre		Post		Pre		Post		Pre		Post	
	1	2	3	4	5	6	7	8	9	10	11	12
A	1:4		1:4		1:4		1:4		1:4		1:4	
B	1:8		1:8		1:8		1:8		1:8		1:8	
C	1:16		1:16		1:16		1:16		1:16		1:16	
D	1:32		1:32		1:32		1:32		1:32		1:32	
E	1:64		1:64		1:64		1:64		1:64		1:64	
F	1:128		1:128		1:128		1:128		1:128		1:128	
G	1:256		1:256		1:256		1:256		1:256		1:256	
H	PCS 1:16		PCS 1:16 + Mn C- Ps 500 ng/ml		C' control		C' control		PCS 1:16 + bacteria only		Medium control	

PCS: positive control serum

5 **A. Plate Quality Control**

Each plate should have its own controls for assay validity. These controls were set up in Row H (wells H1 to H12):

10 Positive serum was diluted in a separate tube to give a final dilution of 8 times the predetermined endpoint titer for positive control serum. When this serum was added to the wells on the plate, the final concentration was 4 times the predetermined endpoint titer. A minimum of 150  $\mu$ l of diluted positive control was needed for each plate.

1. Positive sera control wells (H1 and H2): These wells contained 25 $\mu$ l of positive control serum at a final concentration of 4 times the pre-determined SBA titer
- 15 2. Serotype specific control (H3 and H4): These wells contained 25 $\mu$ l of positive control serum at a final concentration of 4 times the pre-determined SBA titer and 10 $\mu$ l of serotype-specific polysaccharide at 3  $\mu$ g/ml (final on plate concentration is 500 ng/ml).

3. Complement control (i.e., no serum) wells (H5, H6, H7, H8): 25  $\mu$ l of D-PBS-Glc buffer was added.
4. Serum and bacteria only wells (H9 and H10): 12.5  $\mu$ l of D-PBS-Glc buffer was added.
5. Medium alone control (H11 and H12): 50  $\mu$ l D-PBS-Glc buffer was added to these wells (H11 and H12).

Note: When all components were added, the volume of each reagent in Row H control wells was as follows: (complement and bacteria were added in later steps which are also listed here for reference)

	Positive serum	Specificity control	C' control (bacteria +C')	Serum bacteria (serum control)	Medium control
D-PBS-Glc buffer	-	-	25 $\mu$ l	12.5 $\mu$ l	50 $\mu$ l
Polysaccharide	-	10 $\mu$ l	-	-	-
Diluted serum	25 $\mu$ l	25 $\mu$ l	-	25 $\mu$ l	-
C'	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l	-	-
Bacteria	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l	-

10

Reagents were added in the following order:

i) D-PBS-Glc buffer; ii) polysaccharide; and iii) diluted positive control serum.

After adding the above reagents to the control wells, work should be carried out on the test wells.

- 15 Serial dilutions were made of the test serum.

Then dilutions were made of the bacteria.

Bacteria were diluted to the desired dilution.

After making dilutions, an appropriate volume of bacteria at the final dilution was taken out and mixed with equal volume of complement. For one plate, at least 1.2 ml

- 20 diluted bacteria and 1.2 ml complement is recommended.

25  $\mu$ l of bacteria and complement mixture was added to the test and control wells except the "serum and bacteria only" wells and the "medium control" wells.

12.5  $\mu$ l bacteria at final dilution was added to the "serum and bacteria only" wells. For one plate, 50  $\mu$ l diluted bacteria was needed.

- 25 B. Test Sample Dilutions

Note: Serum samples were generally tested without pre-dilution. A minimum of 25  $\mu$ l was determined to be needed to test a serum sample (in duplicate) for bactericidal activity against a single serotype of meningococci. Per specific experiment needs, sera may be heat-inactivated ( $56 \pm 2$  °C,  $30 \pm 5$  minutes) prior to testing. Repeat cycles of freezing and thawing should be avoided. Serum samples were tested at a starting

- 30 dilution of 1:4 (in plate final dilution) and diluted in a two-fold dilution scheme.



1. A sterile reagent reservoir was filled with D-PBS-Glc buffer and a multi-channel pipette was used to add 25  $\mu$ l of D-PBS-Glc buffer to all wells in columns 1-12 except Row A and Row H. Wells in Row H were used for the C' control, medium control and other controls. (See plate template).
2. 37.5  $\mu$ l of D-PBS-Glc buffer was added to all wells in Row A.
3. 12.5  $\mu$ l of undiluted serum was aliquotted into the appropriate wells in Row A according to the plate set-up template.
4. Two-fold serial dilutions of the samples were made using a multichannel pipetter adjusted to 25  $\mu$ l by first transferring 25 $\mu$ l from wells A1 to A12 in the first row to the corresponding wells in Row B; avoiding air bubbles.
5. The contents of the wells were carefully mixed by pipetting up and down ~5 times.
6. 25  $\mu$ l was transferred from the second to the third row and step 5 was repeated.
7. This process was continued on through Row G. The excess 25  $\mu$ l from Row G was discarded into the waste.
8. All wells now contained 25  $\mu$ l. Note: Row A becomes a 1:4 dilution once the remainder of the reagents are added. Row B = 1:8, Row C = 1:16, Row D = 1:32, Row E = 1:64, Row F = 1:128, and Row G = 1:256. Row H contains controls (See plate template above).

## 20

1. A bacteria frozen stock of the desired strain from the  $-70 \pm 15^\circ\text{C}$  freezer was allowed to thaw at room temperature and used immediately; otherwise, the thawed bacteria need be kept on ice until ready to make the dilution, but no longer than 1 hour.
2. The contents of the bacteria vial were mixed by slowly pipetting about one-half the volume of medium in the vial up and down 3 to 5 times; avoiding the introduction of air into the medium. Note: Do not vortex. An initial dilution was made of 1:10 as well as subsequent dilutions needed to reach the desired final concentration. Dilutions were made in D-PBS-Glc buffer (usually in 2-3 steps, e.g., 1:10, then, 1:100 and, finally, 1:30 to make a 1:3000 dilution) in 5 ml round-bottom tubes. The number of bacteria were adjusted to  $\sim 1,000$  cfu/ $12.5\ \mu\text{l}$ , or  $\sim 50$ - $150$  cfu/ $5\ \mu\text{l}$  aliquot from the total volume in the complement control well (usually  $10^3$  to  $10^4$  dilution from working stock if the working stock was prepared when bacteria grew to a density of OD600 around 0.4 to 0.5).

(Frozen aliquots of *N. meningitidis* were pre-tested to determine the appropriate dilution.) Note: Each plate requires ~1.2 ml of diluted bacteria. Dilution varies according to each frozen strain and lot number. Make sure you have a slight (e.g. at least 0.3 ml) excess of diluted bacteria for the assay.

5 D. Bacteria-Complement Addition

1. Frozen stocks of complement were stored at -70°C. The vials of baby rabbit serum were thawed immediately before use. Note: The complement should be at ~4 °C before addition to plates. Two vials of baby rabbit serum that contains at least 1.2 mL combined volume, per assay plate were thawed. All the vials of baby rabbit serum  
10 were combined in a sterile 14 ml tube. Once the complement was pooled, it was added to the plates within 5 minutes after removing the vials from the ice.
2. The appropriate volume of bacteria at the final dilution was taken out, and the bacteria suspension was mixed with an equal volume of complement in a fresh sterile 14 ml tube. A multichannel pipetter set to 25 µl was used to add the bacteria-  
15 complement mix to each of the 96 wells in a microtiter plate except the "serum and bacteria only" control wells and "medium control" wells (H9, H10, H11 and H12). A sterile reagent reservoir was used to hold the bacterial-complement suspension. Where multiple plates were being run, the bacteria-complement suspension was mixed with a pipette, between additions to successive plates.
- 20 3. 12.5 µl bacteria diluted at the final dilution was added to the "serum and bacteria only" control wells.

E. Serum Bactericidal Reaction

- Plates were covered with their lids and the microtiter plates were incubated on a horizontal rotator shaker at 200 rpm at 37 ± 2 °C, No CO<sub>2</sub> for 60 ± 5 min. This step  
25 allowed the killing of bacteria by sera and complement.

F. Viable Bacteria Enumeration

1. Sterile water (RCM-66) was added to the Millipore™ 96 well HV plate, 100µl per well.
2. Filtration: The HV plate was placed on the Millipore™ MultiScreen™ filtration  
30 system and the pump was turned on. Note: Make sure meter is showing 10 ± 2 in. Hg-vac. The plate contents were filtered through the membrane by pushing down the plate with the lid to the base of the unit until the plate was fully emptied and no liquid remained visible.
3. TSB broth (room temperature) was added to the Millipore™ 96 well HV plate;  
35 100µl per well.

4. The contents in the wells from SBA reaction were carefully mixed (4-6 times) and the desired volume of sample (e.g. 5  $\mu$ L) from each well was transferred into the corresponding wells of the HV plates using a multichannel pipette.
- Note: The same volume (either 5  $\mu$ L) must be added to all wells in the HV plates.
- 5 Note: The transferred volume (5  $\mu$ L) is usually 1 /10 of the total volume in each well, except that for the "specificity control" wells, only 5  $\mu$ L is transferred, but the total volume for those wells is 60  $\mu$ L. This disproportion would only affects "serum and bacteria only" wells, won't make any difference on the results of test wells.
5. Filtration: Step 2 was repeated.
- 10 Note: Filtration can happen very fast (within 5 – 10 seconds). Therefore, monitor the filtration carefully to prevent over-drying of the plate. Note that too much broth remaining in the wells will give high background and poor quality staining. Remove plate from the manifold immediately after most of the broth is drawn through the wells, but avoiding over-drying and over-filtration. Carefully monitor wells during
- 15 filtration to observe that well contents have been vacuum-filtered, as discernable by the disappearance of liquid from the outside edges of the wells. Remove plate immediately after this occurs.
6. After filtration, the base of the HV plate unit was pat dry on a paper towel. The plate was wrapped with plastic film.
- 20 7. The wrapped plate was inserted in a plastic bag (to keep moisture inside of plate) and incubated overnight at  $37 \pm 2$  °C,  $5 \pm 1\%$  CO<sub>2</sub>, for a period of 15-17 hours (avoiding the overgrowth of colonies). The plates were placed in an upright position in the incubator.

Volume in each well at this point:

	Positive serum	Specificity control	C' control (bactera +C')	Serum bacteria c Medium (serum control) control
D-PBS-Glc buffer	-	-	25 $\mu$ l	12.5 $\mu$ l
Polysaccharide	-	10 $\mu$ l	-	-
Diluted serum	25 $\mu$ l	25 $\mu$ l	-	25 $\mu$ l
C'	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l	-
Bactena	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l	12.5 $\mu$ l

25

G. Colony staining and Counting

1. Only fresh-filtered stain was used for the HV plates. 0.01% Coomassie blue was filtered with Whatman No. 5 filter paper to remove indissoluble precipitates in the stain. (approximately 5ml of stain required for each HV plate).
- 30 2. Colony staining: The HV plate was placed on the Millipore™ filtration system without pushing down. (Note: the plate is off the vacuum). The pump was turned on.

**Note:** Make sure meter is showing  $10 \pm 2$  in. Hg-vac. 50 $\mu$ l per well of filtered Coomassie blue was added to the plate using the Matrix Multichannel pipette. The plate was only stained for as long as the Coomassie blue was added to whole plate. Then, the plate was pushed down onto the filtration manifold to empty the contents by vacuum filtration.

- 5 3. Destaining: Once the plate was fully emptied, 100 $\mu$ l per well Coomassie blue diluent was immediately added to each well and the plate was kept under the vacuum until no liquid was visible. This step removed most non-specific staining.
- 10 4. The base of the plate was pat dry on a paper towel or equivalent and the base was separated from the 96 well multiscreen plate.
5. The plate was dried upside down in a Biosafety hood with circulating air, until wells were completely dried (approximately  $10 \pm$  minutes)
6. ImmunoSpot™ counter setting for colony counting

	Standard setting	
	(medium size colonies)	(small colonies)
15 Sensitivity:	100	140
Spot size:	0.12	0.06
Separation tolerance:	10.0	10.0
Diffuse counting is:	on	on
20 Diffuseness:	20	20
Overdeveloped area handling is:	off	off
Color system:	Blue	Blue
Detailed counting is:	off	off
Background balance is:	on	on
25 Background balance:	110	40
Hole filling is:	on	off

**Note:** ImmunoSpot™ settings should be calibrated and adjusted (at 3 month intervals) by comparison with counts obtained by human readers.

The longer the overnight incubation, the bigger the spot size.

- 30 7. The plate was scanned in the imaging system, the image was printed and the data was counted. The image data were saved as tiff files on CD-R disk.

#### H. SBA titer calculations

The bactericidal titer is the reciprocal of the serum dilution with at least 50% killing, when compared to the average growth in the complement control wells. Occasionally,

sera with high titers need to be re-tested at higher initial dilutions than 1:4 (for example, starting from 1:64) to determine the phagocytic titer. Serum samples with phagocytic titers <4 are reported as a titer of 2 for purposes of data analysis.

#### *I. Experimental Results*

- 5 SBA assays were individually run with bacterial samples of *N. meningitidis* serotypes A, B, C, Y, and W135, and *E. coli*. Bacterial growth or lack thereof (due to the effect of an antibiotic substance) was able to be detected and analyzed within the multi-well plates. Figure 4 shows the basic assay template set-up utilized to run the experiment with serum samples of *N. meningitidis* serotype C. Figures 5 and 6, respectively,  
10 show the colony read-out obtained from the wells both pictorially and numerically.

#### EXAMPLE 5. PROCEDURE FOR VALIDATION OF OPSONOPHAGOCYTTIC ASSAY

- Analytical validation of the OPK assay for *Streptococcus pneumoniae* serotype 23F  
15 was performed in six assay runs utilizing two different operators (LS and SW), and three cell batches (p38, p39, and p71). Each of the six runs was performed using two types of plates (agar and HV). The agar plates were manually counted whereas the HV plates were counted both manually and via automation under "standard" and "high" sensitivity. Test samples included in the validation consisted of 48 ELISA  
20 negative samples, and three pools of pediatric sera ranging from low/negative to high OPK response as assessed in preliminary runs. The serum pools were tested within each run while the negative samples were evenly divided across runs. Each plate also included (1) four "No Serum" control wells containing bacteria, complement, and HL-60 cells but no antisera; (2) two "Medium Only" control wells, (3) a positive control  
25 sera (QC-1) tested at three dilutions, and (4) two specificity controls (QC-1 serum tested with 23F PS and QC-1 serum tested with C-PS). The experimental layout is provided in Figure 7. Specificity was also assessed by determining the ability of polysaccharides of a known serotype (6B, 9V, 14, 18C, 19F, 23F, and C-Ps at 1 µg/ml) to inhibit killing by positive control antisera (Pool 019 and QC-1 tested at the  
30 1:64 dilution). The experimental layout for assessing specificity is provided in Figure 8. OPK titers were determined by running serial 2-fold dilution of serum in duplicate assay wells.

The results of the analytical validation are summarized below:

Assay Characteristic	Validation Result
Count Variability	The automated counting method (standard sensitivity) was more consistent than the manual method. Within the HV plate, the inter-run and intra-run %RSD's were 1.7% and 2.9%, respectively, for automated counting as compared to 14.3% and 14.2%, respectively, for manual counting.
Relationship Between Count Level and Volume Plated	The relationship between colony count and volume is linear between 5µl and 15µl. For every 2-fold increase in volume, the number of colonies is increased by 1.59-fold (95% CI = (1.33, 1.89)).
CFU Level Comparison: Periphery Wells versus Central Wells	Relative to CFU levels in the central wells, CFU levels in the periphery wells were reduced by 3.1% on average (95% CI=(0.2, 9 CFU)). The percent reduction in CFU levels tended to increase with increasing plate volume as the reductions at 5ul, 10ul, and 15ul were 1.1% , 3.3%, and 4.3%, respectively.
Extra-variability	Extra-variability is assessed on the square roots of the replicate counts. The specification limit, applied to the range of the transformed counts, is 2.57.
Ruggedness (Cell Passage, Operator, Plate, Counting Method)	Differences in OPK titer between cell passages (p38-p71), operators, plates and counting methods were all within 2-fold on average, and therefore the assay is judged to be acceptably rugged to changes in these parameters. There was no evidence that the counting method within the HV plate (Manual, Standard, and High) affected titer. While the HV and Agar plates resulted in similar titer levels on average, individual results tended to be more variable on the Agar plate than on the HV plate.

Assay Characteristic	Validation Result
Precision	Variability for the HV plates did not exceed that for the Agar plates. With respect to reported titers, there was no evidence of a difference in assay precision between the manual and automated counting methods within the HV plate. The assay precision (percent relative standard deviation in OPK titer) for the HV plates is approximately 70% which is acceptable for an endpoint dilution assay.
Sero-Status Cutoff	The titer distribution of the negative control samples conservatively supports a serostatus cutoff of 1:8 in the OPK assay. Thus, samples having a titer of <1:8 are considered negative in the assay.
Specificity	Of the seven polysaccharides evaluated (6B, 9V, 14, 18C, 19F, 23F, and C-Ps at 1 µg/ml), only 23F inhibited killing.
Assay Controls	The following six requirements are recommended in order to qualify an assay run: (1) the positive control must be positive (i.e., >50% inhibition) at least at one of three dilutions tested, 1:0.5X, 1:X, and 1:2X, where X represents the endpoint titer for the positive control; (2) The positive control tested at a dilution of 1:0.5X in the presence of 23F-Ps at 1 µg/ml must be negative (i.e., ≤50% inhibition); (3) The positive control tested at a dilution of 1:0.5X in the presence of C-Ps at 1 µg/ml must be positive (i.e., >50% inhibition); (4) The replicate counts within the "Medium Only" and "No Serum" control wells must satisfy the extra-variability criteria; (5) The average number of colonies across the four "No Serum" control wells must be ≥25 after correcting for the "Medium Only" control wells and <200 (<300) for Agar (HV) plates prior to correcting for the "Medium Only" control wells; and (6) The average number of colonies in the two "Medium Only" control wells must be ≤25% of the average number of CFU in the four "No Serum" control wells. Any exceptions would require

	appropriate notation in order to accept the run.

**EXAMPLE 6. SBA ASSAY APPLICATION: Rapid Selection of Functional Monoclonal Antibodies that kill *Neisseria Meningitidis* Serotype B**

- The selection of hybridomas that produce monoclonal antibodies typically involves the screening of hundreds of culture supernatants for antibodies of the desired specificity and function. Usually, this involves a two-step selection process in which hybridomas are screened for the ability to produce antibodies that bind to an antigen in an enzyme-based immunoassay (EIA). Subsequently, the supernatants are screened for functional activity (e.g. utilizing an OPK or SBA assay). Utilizing the methods described in the instant application, these separate steps are combined into a single rapid screening assay. Results employing such an assay for the screening of hybridoma culture supernatants for serum bactericidal activity against *N. meningitidis* serotype B are exemplified in Figure 9. In the Figure, one can clearly see serum bactericidal activity in the hybridoma culture supernatants of wells F4 and F11 which, incidentally, have only a few bacterial colonies as opposed to the remaining wells that all have a large numbers of colonies. It is important to note that antibiotics must not be added to hybridoma culture medium as they will interfere with the rapid SBA screening procedure by killing the test bacteria.



## WHAT IS CLAIMED IS:

1. A method for enumerating microbial colonies in a sample comprising:
  - (a) transferring a sample comprising a microbe(s) of interest in a liquid  
5 medium to the wells of a multi-well filter plate;
  - (b) removing excess media from the wells;
  - (c) allowing sufficient time for the microbe(s) to grow into discrete colonies on residual growth media captured within and under the filter plate; and
  - (d) enumerating the microbial colonies in the sample by using a means  
10 suitable for enumeration of samples in multi-well format.
2. A method in accordance with claim 1 wherein the microbe(s) are bacteria and the microbial colonies are bacterial colonies.
3. A method in accordance with claim 1 wherein the microbe(s) are yeast and the microbial colonies are colonies of yeast.
- 15 4. A method in accordance with claim 1 wherein the microbe(s) are fungi and the microbial colonies are fungal colonies.
5. A method in accordance with claim 1 wherein the multi-well filter plate of step (a) comprises growth medium.
6. A method in accordance with claim 1 wherein the multi-well filter plate is a 96  
20 well filter plate.
7. A method in accordance with claim 2 wherein the bacteria are grown on the filter plate for a period of 14-18 hours.
8. A method in accordance with claim 1 wherein the filter plate is a Millipore™ 96 well HV plate.
- 25 9. A method in accordance with claim 1 wherein the filter plate is a Millipore™ Multiscreen™ HV 0.45 µm Opaque Sterile Filtration plate.
10. A method in accordance with claim 1 wherein the excess media is removed by vacuum filtration.
11. A method in accordance with claim 1 wherein the excess media is removed by  
30 centrifugation.
12. A method in accordance with claim 1 wherein the enumeration of microbial colonies is accomplished with a device capable of acquiring images and/or information from wells in multi-well format.
13. A method in accordance with claim 12 wherein the device is capable of  
35 acquiring images and/or information from wells in 96 well format.

14. A method in accordance with claim 12 wherein the number of bacteria in the sample is determined using a computer-assisted video imaging and analysis system
15. A method in accordance with claim 12 wherein the number of bacteria in the sample is determined using an ImmunoSpot™ Analyzer.
- 5 16. A method in accordance with claim 1 wherein the microbe(s) prior to step (a) had been contacted with an antimicrobial agent.
17. A method in accordance with claim 16 wherein the microbe(s) prior to step (a) had been contacted with an antimicrobial antiserum.
18. A method in accordance with claim 16 wherein the microbe(s) prior to step (a) were further contacted with complement or active components thereof.
- 10 19. A method in accordance with claim 17 wherein the microbe(s) prior to step (a) were further contacted with an effector cell capable of engulfing the microbe(s).
20. A method in accordance with claim 18 wherein the effector cell is a phagocyte.
- 15 21. A method in accordance with claim 18 wherein the effector cell is a differentiated HL-60 cell.
22. A method in accordance with claim 18 wherein the effector cell is a peripheral blood polymorphonuclear leukocyte.
23. A method in accordance with claim 2 wherein the bacteria is a gram-positive bacteria.
- 20 24. A method in accordance with claim 2 wherein the bacteria is a gram-negative bacteria.
25. A method in accordance with claim 2 wherein the bacteria is a pathogenic microorganism.
- 25 26. A method in accordance with claim 2 wherein the bacteria is *Streptococcus pneumoniae*.
27. A method in accordance with claim 2 wherein the bacteria is *Neisseria meningitidis*.
28. A method in accordance with claim 2 wherein the bacteria is *Escherichia coli*.
- 30 29. A method in accordance with claim 2 wherein the bacteria is *Staphylococcus aureus*.
30. A method in accordance with claim 2 wherein the bacteria is *Bacillus anthracis*.
31. A method for analyzing microbe(s), their growth and/or viability in a sample comprising:
- 35

32. (a) transferring a sample comprising a microbe(s) of interest in a liquid medium to the wells of a multi-well filter plate;
- (b) removing excess media from the wells;
- 5 (c) allowing sufficient time for the microbe(s) to grow into discrete colonies on residual growth media captured within and under the filter plate; and
- (d) analyzing the microbe(s), their growth and/or viability in the sample by a means suitable for analysis of samples in multi-well format.
33. A method in accordance with claim 30 wherein the microbe(s) are bacteria.
- 10 34. A method in accordance with claim 30 wherein the microbe(s) are yeast.
35. A method in accordance with claim 30 wherein the microbe(s) are fungi.
36. A method in accordance with claim 30 wherein the multi-well filter plate of step (a) comprises growth medium.
37. A method in accordance with claim 34 wherein the multi-well filter plate is a 15 96 well filter plate.
38. A method in accordance with claim 31 wherein the bacteria are grown for a period of 14-18 hours.
39. A method in accordance with claim 30 wherein the filter plate is a Millipore™ 96 well HV plate.
- 20 40. A method in accordance with claim 30 wherein the filter plate is a Millipore™ Multiscreen™ HV 0.45 µm Opaque Sterile Filtration plate.
41. A method in accordance with claim 30 wherein the excess media is removed by vacuum filtration.
42. A method in accordance with claim 30 wherein the excess media is removed 25 by centrifugation.
43. A method in accordance with claim 30 wherein microbe(s), their growth and/or viability is analyzed with a device capable of acquiring images and/or information from wells in multi-well format.
44. A method in accordance with claim 41 wherein the device is capable of 30 acquiring images and/or information from wells in 96 well format.
45. A method in accordance with claim 42 wherein the number of bacteria in the sample is determined using a computer-assisted video imaging and analysis system.
46. A method in accordance with claim 42 wherein the microbe(s) are analyzed 35 using an ImmunoSpot™ Analyzer.

47. A method in accordance with claim 30 wherein the microbe(s) prior to step (a) had been contacted with an antimicrobial agent.
48. A method in accordance with claim 30 wherein the microbe(s) prior to step (a) had been contacted with antimicrobial antiserum.
- 5 49. A method in accordance with claim 45 wherein the microbe(s) prior to step (a) were further contacted with complement or active components thereof.
50. A method in accordance with claim 46 wherein the microbe(s) prior to step (a) were further contacted with an effector cell capable of engulfing the microbe(s).
- 10 51. A method in accordance with claim 47 wherein the effector cell is a phagocyte.
52. A method in accordance with claim 47 wherein the effector cell is a differentiated HL-60 cell.
53. A method in accordance with claim 47 wherein the effector cell is a peripheral blood polymorphonuclear leukocyte.
- 15 54. A method in accordance with claim 31 wherein the bacteria is a gram-positive bacteria.
55. A method in accordance with claim 31 wherein the bacteria is a gram-negative bacteria.
56. A method in accordance with claim 31 wherein the bacteria is a pathogenic microorganism.
- 20 57. A method in accordance with claim 31 wherein the bacteria is *Streptococcus pneumoniae*.
58. A method in accordance with claim 31 wherein the bacteria is *Neisseria meningitidis*.
- 25 59. A method in accordance with claim 31 wherein the bacteria is *Escherichia coli*.
60. A method in accordance with claim 31 wherein the bacteria is *Staphylococcus aureus*.
61. A method in accordance with claim 31 wherein the bacteria is *Bacillus anthracis*.
- 30 62. A method for evaluating antimicrobial agents comprising:
- (a) contacting a sample comprising a microbe(s) of interest with the antimicrobial agent;
- (a) transferring the sample comprising the microbe(s) of interest in a liquid medium to the wells of a multi-well filter plate;
- 35

- (b) removing excess media from the wells;
  - (c) allowing sufficient time for the microbe(s) to grow into discrete colonies on residual growth media captured within and under the filter plate; and
  - 5 (d) evaluating the effect of the antimicrobial agent on the growth and/or viability of the microbe(s) with a means suitable for analysis of samples in multi-well format.
63. A method for evaluating antimicrobial agents comprising:
- 10 (a) transferring a sample comprising a microbe(s) of interest in a liquid medium to the wells of a multi-well filter plate;
  - (b) contacting the sample comprising the microbe(s) of interest with the antimicrobial agent;
  - (c) removing excess media from the wells;
  - (d) allowing sufficient time for the microbe(s) to grow into discrete colonies on residual growth media captured within and under the filter plate; and
  - 15 (e) evaluating the effect of the antimicrobial agent on the growth and/or viability of the microbe(s) with a means suitable for analysis of samples in multi-well format.
- 20 64. A method in accordance with claim 60 wherein the antimicrobial agent is a monoclonal antibodies present within hybridoma culture supernatant.

## TITLE OF THE INVENTION

NOVEL PROCEDURE FOR GROWTH, IMAGING, AND ENUMERATION OF  
MICROBIAL COLONIES FOR SEROLOGICAL OR SCREENING ASSAYS

## 5 ABSTRACT OF THE DISCLOSURE

A novel method is provided, for use within serological or screening assays, wherein microbial colonies are grown on filter membranes in multi-well plates. This process enables the colonies to be stained, imaged, and counted automatically using such automated systems as, *e.g.*, computer and video-based

10 imaging systems.

**FIG. 1**

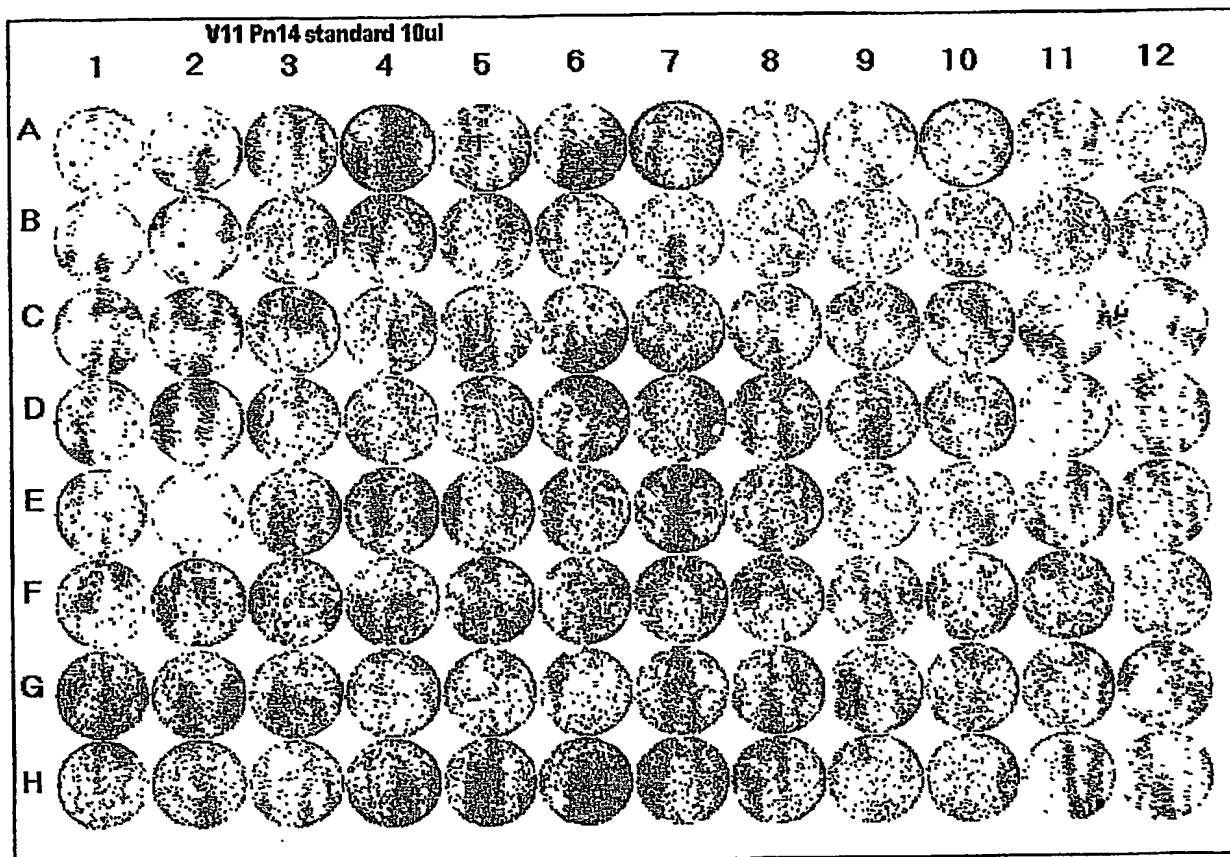


FIG. 2



DPK V11-14 Colony Read Out by Immunospot Counter												
	(019 pool)		(014 pool)		(013 pool)		Neg. Sera & Dilution		Neg. Sera & Dilution		Controls	
	21	16	61	51	29	33	185	165	183	203	133	128
	8	15	74	118	95	107	145	194	178	168	137	160
	15	9	97	93	101	86	174	197	170	198	8	14
	19	9	95	107	106	98	165	186	158	200	17	14
	19	7	101	109	105	120	166	177	94	101	25	42
	68	42	95	104	142	105	165	179	123	98	161	187
	136	145	121	125	92	105	183	198	168	158	43	21
	147	133	128	127	114	109	178	168	182	180	9	5

**FIG. 3**

Serotype C					Serotype C bacteria, MJC post serum 1:16						
Pre	Post			Controls			Post		Post		
1	2	3	4	5	6	7	8	9	10	11	12
A 1:4		1:4		Cells only (2) pre 1:16		Mn C-Ps 10000 ng/ml	W135-Ps 10000 ng/ml	Y-Ps 10000 ng/ml			
B 1:8		1:8		Cells only (2) post 1:16		Mn C-Ps 1000 ng/ml	W135-Ps 1000 ng/ml	Y-Ps 1000 ng/ml			
C 1:16		1:16		Cells + C'		Mn C-Ps 100 ng/ml	W135-Ps 100 ng/ml	Y-Ps 100 ng/ml			
D 1:32		1:32				Mn C-Ps 10 ng/ml	W135-Ps 10 ng/ml	Y-Ps 10 ng/ml			
E 1:64		1:64		post 1:16 + Mn C-PS 10 ug/ml (2)		post 1:16 no Ps	post 1:16 no Ps	post 1:16 no Ps			
F 1:128		1:128		post 1:16 + Mn W135-PS 10 ug/ml		Cells only post 1:16	Cells only post 1:16	Cells only post 1:16			
G 1:256		1:256		post 1:16 + Mn Y-PS 10 ug/ml (2)		Cells+C'	Cells+C'	Cells+C'			
H 1:512		1:512		medium control		medium control	medium control	medium control			

FIG. 4

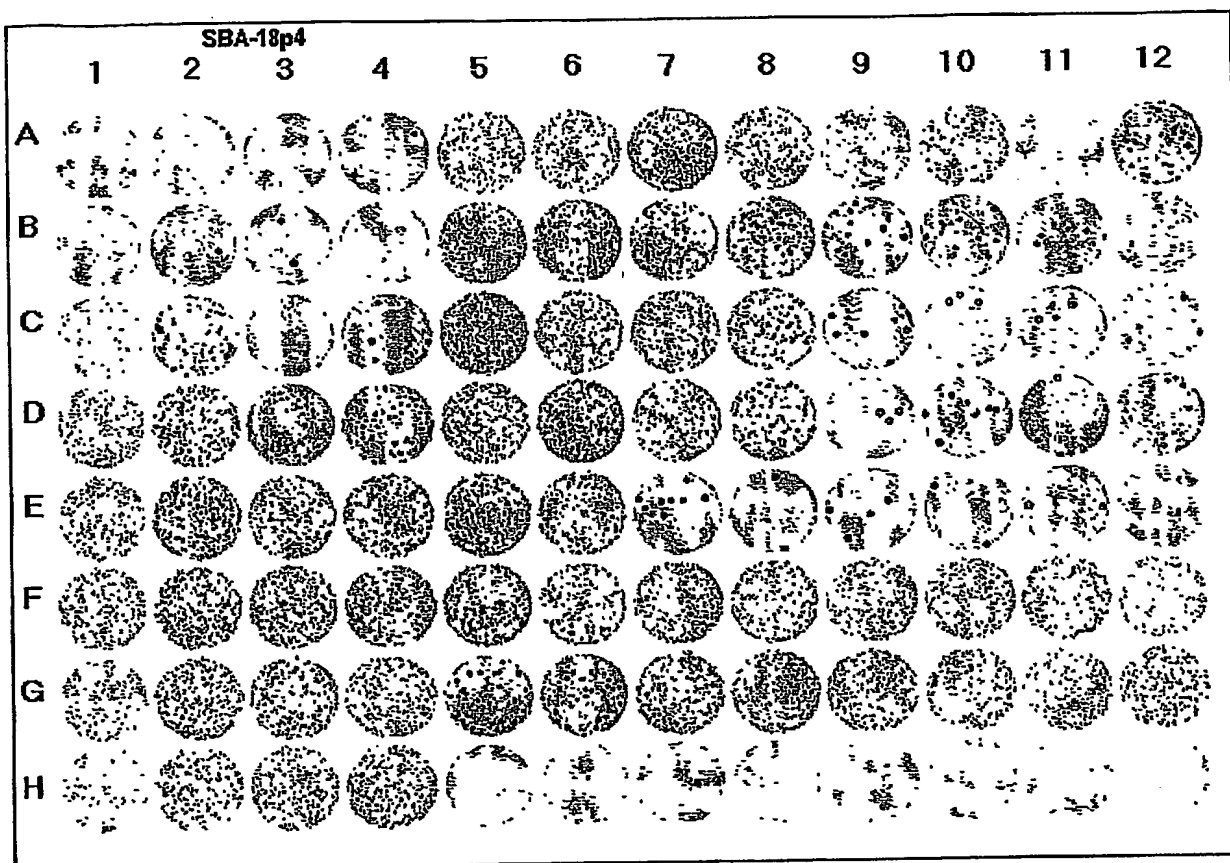


FIG. 5

Serotype C										Serotype C bacteria, MJC post serum 1:16									
pre			post		controls			Mn C-ps			Mn W135-ps			Mn Y-ps					
	Avg	%Killing		Avg	%Killing		Avg	%Killing		Avg	%Killing		Avg	%Killing		Avg	%Killing		
A	6	96	1	100	3	152	7	Cells only pre 1:16	Mn C-Ps 10000 ng/ml	121	7	W135-Ps 10000 ng/ml	73	47	Y-Ps 10000 ng/ml	34	78		
B	6	96	1	99	16	131	26	Cells only post 1:16	Mn C-Ps 1000 ng/ml	96	26	W135-Ps 1000 ng/ml	29	79	Y-Ps 1000 ng/ml	9	94		
C	50	68	11	93		156	14	Cells + C (4 wells)	Mn C-Ps 100 ng/ml	112	14	W135-Ps 100 ng/ml	7	95	Y-Ps 100 ng/ml	9	94		
D	131	16	70	55			24		Mn C-Ps 10 ng/ml	99	24	W135-Ps 10 ng/ml	10	93	Y-Ps 10 ng/ml	29	81		
E	126	19	112	29	32	107	91	post 1:16 Mn C-Ps (2)	post 1:16 no Ps	12	91	post 1:16 no Ps	6	96	post 1:16 no Ps	8	95		
F	146	7	128	18	38	97	152	post 1:16 Mn W135-Ps 10 ug/ml (2)	Cells only post 1:16	152	-17	Cells only post 1:16	155	-14	Cells only post 1:16	132	12		
G	167	-7	129	18	71	45	130	post 1:16 Mn Y-Ps 10 ug/ml (2)	Cells + C	130	0	Cells + C	136	0	Cells + C	150	0		
H	173	-11	161	-3	100	0	0	medium control	medium control	0	100	medium control	0	100	medium control	0	100		

FIG. 6

**Attachment I**  
**Plate Layout for Each Validation Run**  
**Each Run was Performed Using the Agar and HV Plate Methods**

Serotype 23F

Plate Operator HL-60  
 1 1 1 Batch

(019 pool)		(014 pool)		(013 pool)		neg. sera		neg. sera		Controls	
1	2	1	2	1	2	#	Dilution	#	Dilution	1	2
1:32		1:8		1:8		1	1:8	5	1:8	Cells + C'	
1:64		1:16		1:16			1:16		1:16	(no serum)	
1:128		1:32		1:32		2	1:8	6	1:8	QC-1 (2X titer)	
1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
1:512		1:128		1:128		3	1:8	7	1:8	QC-1 (0.5X titer)	
1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
1:2048		1:512		1:512		4	1:8	8	1:8	QC-1 (1:64) + C-Ps*	
1:4096		1:1024		1:1024			1:16		1:16	Medium alone	

Note: "2X titer" = 50% endpoint titer of positive control serum times 2.  
 \*Use polysaccharides (PS) at 1 mcg/mL

Plate Operator HL-60  
 2 2 1 Batch

(019 pool)		(014 pool)		013 pool		neg. sera		neg. sera		Controls	
1	2	1	2	1	2	#	Dilution	#	Dilution	1	2
1:32		1:8		1:8		9	1:8	13	1:8	Cells + C'	
1:64		1:16		1:16			1:16		1:16	(no serum)	
1:128		1:32		1:32		10	1:8	14	1:8	QC-1 (2X titer)	
1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
1:512		1:128		1:128		11	1:8	15	1:8	QC-1 (0.5X titer)	
1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
1:2048		1:512		1:512		12	1:8	16	1:8	QC-1 (1:64) + C-Ps*	
1:4096		1:1024		1:1024			1:16		1:16	Medium alone	

Plate Operator HL-60  
 3 1 2 Batch

(019 pool)		(014 pool)		013 pool		neg. sera		neg. sera		Controls	
1	2	1	2	1	2	#	Dilution	#	Dilution	1	2
1:32		1:8		1:8		17	1:8	21	1:8	Cells + C'	
1:64		1:16		1:16			1:16		1:16	(no serum)	
1:128		1:32		1:32		18	1:8	22	1:8	QC-1 (2X titer)	
1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
1:512		1:128		1:128		19	1:8	23	1:8	QC-1 (0.5X titer)	
1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
1:2048		1:512		1:512		20	1:8	24	1:8	QC-1 (1:64) + C-Ps*	
1:4096		1:1024		1:1024			1:16		1:16	Medium alone	

FIG. 7A

**Attachment I (continued)**  
**Plate Layout for Each Validation Run**  
**Each Run was Performed Using the Agar and HV Plate Methods**

Plate	Operator	HL-60 Batch	(019 pool)		(014 pool)		013 pool		neg. sera		neg. sera		Controls	
			1	2	1	2	1	2	#	Dilution	#	Dilution	1	2
4	2	2	1:32		1:8		1:8		25	1:8	29	1:8	Cells + C'	
			1:64		1:16		1:16			1:16		1:16	(no serum)	
			1:128		1:32		1:32		26	1:8	30	1:8	QC-1 (2X titer)	
			1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
			1:512		1:128		1:128		27	1:8	31	1:8	QC-1 (0.5X titer)	
			1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
			1:2048		1:512		1:512		28	1:8	32	1:8	QC-1 (1:64) + C-Ps*	
			1:4096		1:1024		1:1024			1:16		1:16	Medium alone	
5	1	3	1:32		1:8		1:8		33	1:8	37	1:8	Cells + C'	
			1:64		1:16		1:16			1:16		1:16	(no serum)	
			1:128		1:32		1:32		34	1:8	38	1:8	QC-1 (2X titer)	
			1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
			1:512		1:128		1:128		35	1:8	39	1:8	QC-1 (0.5X titer)	
			1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
			1:2048		1:512		1:512		36	1:8	40	1:8	QC-1 (1:64) + C-Ps*	
			1:4096		1:1024		1:1024			1:16		1:16	Medium alone	
6	2	3	1:32		1:8		1:8		41	1:8	45	1:8	Cells + C'	
			1:64		1:16		1:16			1:16		1:16	(no serum)	
			1:128		1:32		1:32		42	1:8	46	1:8	QC-1 (2X titer)	
			1:256		1:64		1:64			1:16		1:16	QC-1 (1X titer)	
			1:512		1:128		1:128		43	1:8	47	1:8	QC-1 (0.5X titer)	
			1:1024		1:256		1:256			1:16		1:16	QC-1 (1:64) + 23F PS	
			1:2048		1:512		1:512		44	1:8	48	1:8	QC-1 (1:64) + C-Ps*	
			1:4096		1:1024		1:1024			1:16		1:16	Medium alone	

FIG. 7B

**Attachment II**  
**Plate Layout for Specificity Experiment**  
 One Agar Plate and One HV Plate was Prepared by Operator SW

Column:	1	2	3	4	5	6	7	8	9	10	11	12
Serum:	Pool 019		QC-1		Pool 019 minus C'		QC-1 minus C'		Pool 019 Specificity		QC-1 Specificity	
Row	1	2	1	2	1	2	1	2	1	2	1	2
A	1:32		1:32		1:32		1:32		1:64		1:64	
B	1:64		1:64		1:64		1:64		1:64 + Ps 6B		1:64 + Ps 6B	
C	1:128		1:128		1:128		1:128		1:64 + Ps 9V		1:64 + Ps 9V	
D	1:256		1:256		1:256		1:256		1:64 + Ps 14		1:64 + Ps 14	
E	1:512		1:512		1:512		1:512		1:64 + Ps 18C		1:64 + Ps 18C	
F	1:1024		1:1024		1:1024		1:1024		1:64 + Ps 19F		1:64 + Ps 19F	
G	1:2048		1:2048		1:2048		1:2048		1:64 + Ps 23F		1:64 + Ps 23F	
H	No Serum		No Serum		No Serum		No Serum		1:64 + C-Ps		1:64 + C-Ps	

FIG. 8

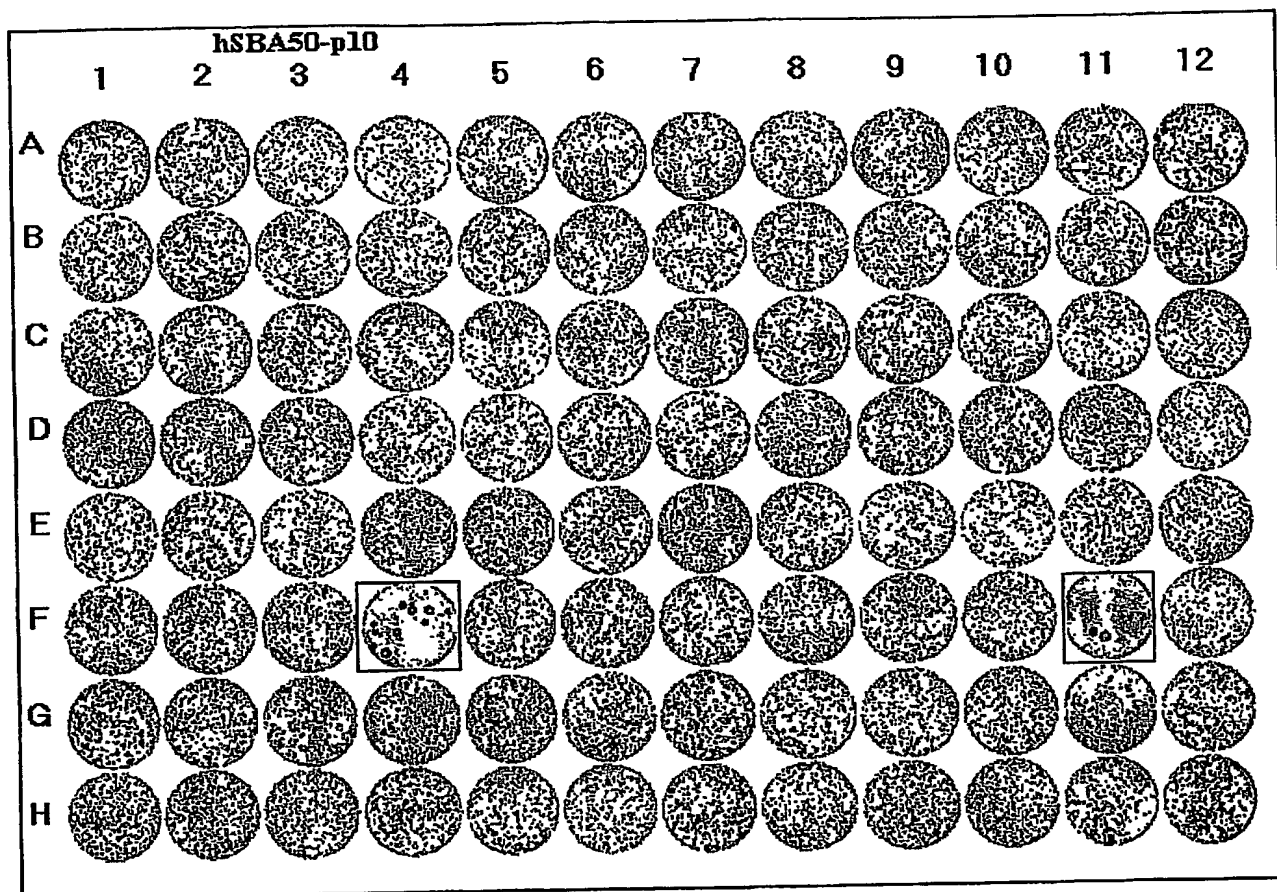


FIG. 9



**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**